


Containment Drum-Dispensing into drums or small bags with containment level OEB 5

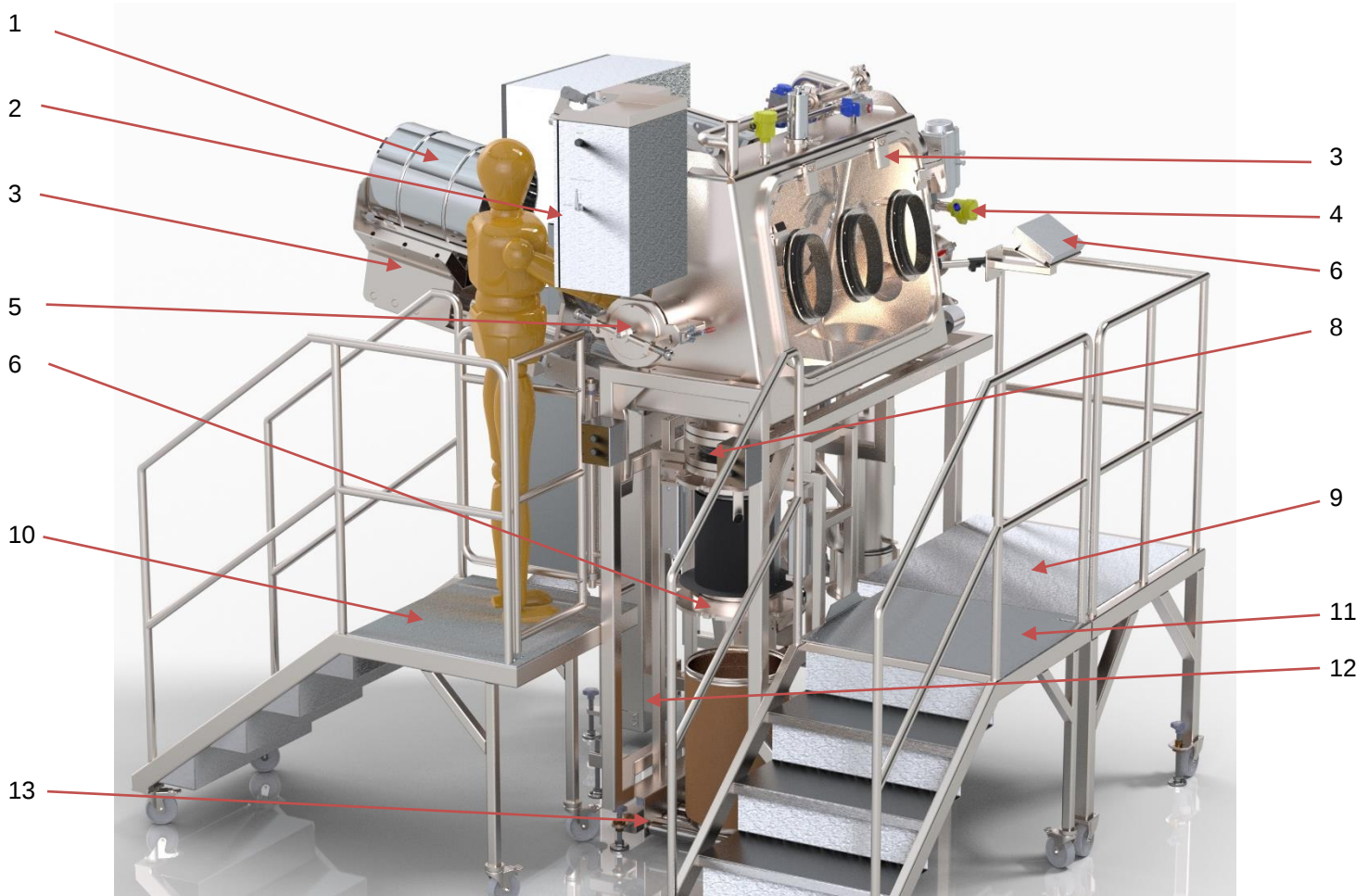


1. Drum in docked position to the isolator. Docking with SmartDock
2. Operating panel for docking / undocking drum liner
3. Pneumatic drive movement to move the tilted drum to, and away from the SmartDock
4. Isolator with hinged glass door
5. Bag out port to remove the shower cap of the last drum
6. Drum filling with continuous liner.
7. Weighing panel of scale in the isolator
8. Decoupling compensator between Isolator and BigBag-filling system which is on a scale
9. Platform part guided on a rail to access bag docking for small amounts
10. Mobile platform for drum docking / undocking
11. Mobile platform
12. Stationary lifting device with automatic drives for
 - 12.1. Lifting / Lowering
 - 12.2. Step less tilting
 - 12.3. Horizontal movement of the drum
13. Balance

- Containment Docking
- Solids calibration
- Lift-and Tilting devices
- Container-and Drum Blending

CH-4431 Bennwil
Tel. +41-(0)61-956 91 91
Fax +41-(0)61-956 91 90
info@rubitec.ch
www.rubitec.ch

rubitec ag

Containment Drum-Dispensing into drums or small bags with containment level OEB 5

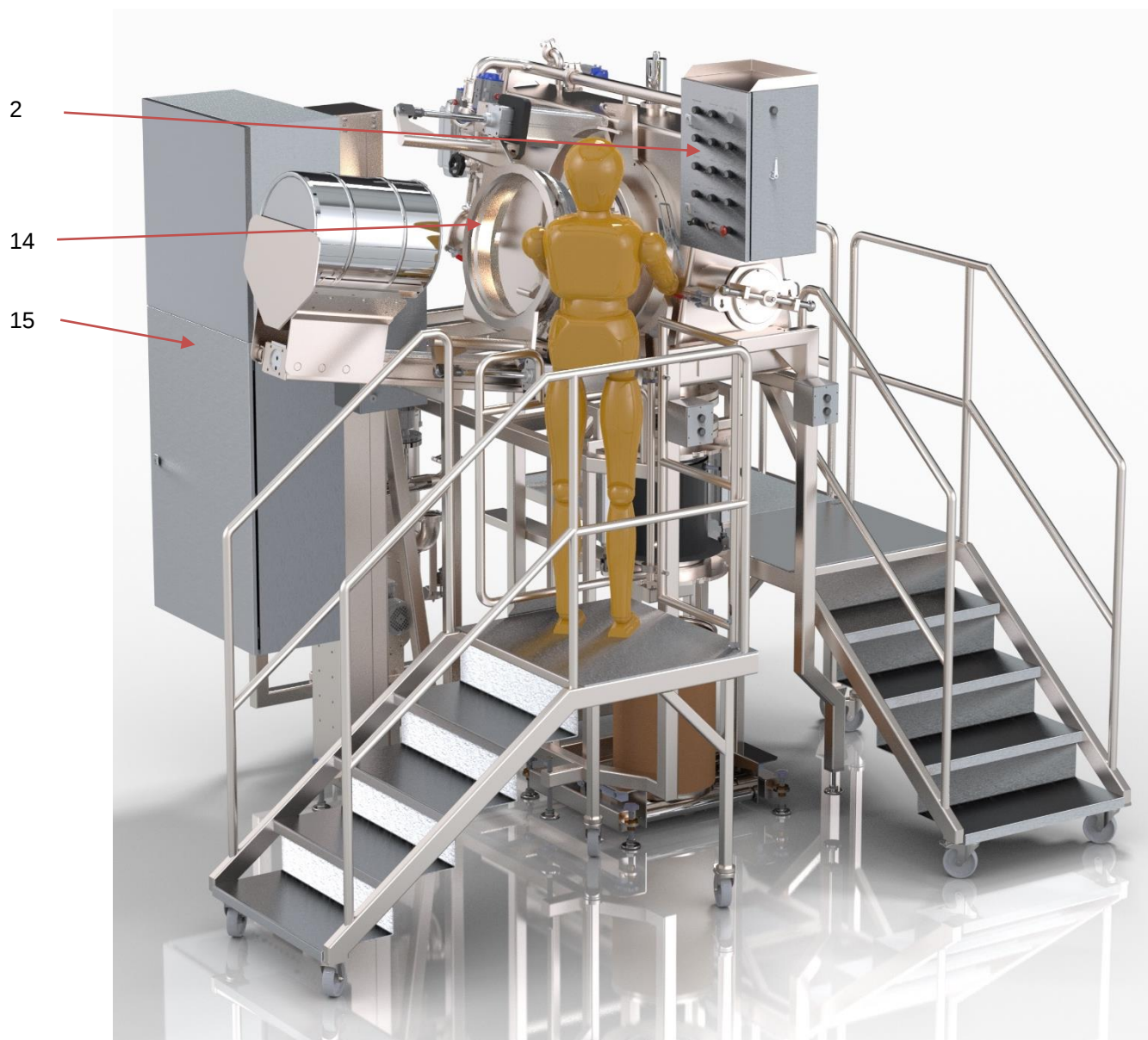

- 2 Operating panel for docking / undocking drum liner
- 14 SmartDock (patented liner docking system)
- 15 Control cabinet

- Containment Docking
- Solids calibration
- Lift-and Tilting devices
- Container-and Drum Blending

CH-4431 Benwil
Tel. +41-(0)61-956 91 91
Fax +41-(0)61-956 91 90
info@rubitec.ch
www.rubitec.ch

rubitec ag